

NAME: .....SIGNATURE: .....

P525/3

CHEMISTRY

(Practical)

Paper 3

Aug/Sept. 2024

3  $\frac{1}{4}$  hours.

## INTERNAL MOCK EXAMS

2024

S6 CHEMISTRY PRACTICAL

Paper 3

3 hours 15 minutes

### INSTRUCTIONS:

Answer **all** questions. Use **blue** or **black** ink. Any work done in pencil will not be marked **except** drawings.

Record your answers on this question paper in the spaces provided.

Mathematical tables and silent non-programmable scientific calculators may be used.

Reference books (i.e textbooks, booklets on qualitative analysis, etc) should **not** be used.

Candidates are not allowed to start working with the apparatus for the first **15 minutes**. This time is to enable candidates read the question paper and make sure they have all the apparatus and chemicals that they may need.

| For Examiners' Use only |     |     |       |
|-------------------------|-----|-----|-------|
| Q.1                     | Q.2 | Q.3 | Total |
|                         |     |     |       |
|                         |     |     |       |

1. You are provided with the following:

**FA1**, which is an aqueous solution of a strong acid,  $H_2X$ .

**FA2**, which is a solution containing  $4.25 \text{ g l}^{-1}$  of hydroxide ions,  $\bar{O}H$ .

**FA3**, which is a  $0.04M$  solution of hydroxide ions.


Solid **Z** which is a basic compound.

**You are required to determine the;**

- (i) concentration of **FA1** in moles per litre
- (ii) stoichiometric ratio of reaction between  $H_2X$  in **FA1** and **Z**.

**Theory;**

Hydroxide ions react with hydrogen ions according to the following equation:

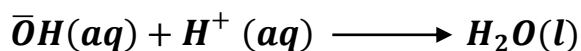

**Procedure A;**

- (a) Using a measuring cylinder, transfer  $40 \text{ cm}^3$  of **FA1** into a  $250 \text{ cm}^3$  volumetric flask and make up to the mark with distilled water. Label the solution formed **FA4**.
- (b) Pipette  $20.0$  (or  $25.0$ )  $\text{cm}^3$  of **FA2** into a clean conical flask, add 2-3 drops of phenolphthaein indicator and titrate it with **FA4** from the burette.
- (c) Repeat the titration until you obtain consistent results.
- (d) Record your results in **Table 1** below.

**Results:**

**Table 1**

Volume of pipette used;.....  $\text{cm}^3$  (  $\frac{1}{2}$  mark)

|  |  |  |  |
|--|--|--|--|
| Final burette reading( $\text{cm}^3$ )     |  |  |  |
| Initial burette reading( $\text{cm}^3$ )   |  |  |  |
| Volume of <b>FA4</b> used( $\text{cm}^3$ ) |  |  |  |

(4  $\frac{1}{2}$  marks)

Titre values used for calculating average volume of **FA4**.

( $\frac{1}{2}$  marks)

.....cm<sup>3</sup>

used.....  $\text{cm}^3$  (2  $\frac{1}{2}$  marks)

Calculate the;

(i) number of moles of  $H_2X$  in  $250\text{cm}^3$  of **FA4**. (4  $\frac{1}{2}$  marks)  
( $H = 1$ ;  $O = 16$ )

[illegible]

(ii) concentration of  $H_2X$  in moles per litre of **FA1**. (1  $\frac{1}{2}$  marks)

.....

.....

.....

.....

.....

**Procedure B;**

- d) Weigh accurately about 2.5 g of **Z** into a clean beaker. Using a measuring cylinder, add to it 20.0cm<sup>3</sup> of **FA1** and stir to dissolve. Transfer the contents of the beaker into a 250 cm<sup>3</sup> volumetric flask and add more distilled water to top up to the mark. Label the resultant solution **FA5**.
- e) Pipette 20.0 (or 25.0)cm<sup>3</sup> of **FA5** into a clean conical flask, add 2-3 drops of phenolphthalein indicator and titrate it with **FA3** from the burette.
- f) Repeat the titration until you obtain consistent results.
- g) Record your results in **Table 2** below.

**Results:**

( 1  $\frac{1}{2}$  marks)

Mass of **Z** and weighing bottle:.....g

Mass of empty weighing bottle:.....g

Mass of **Z** used:.....g

**Table 2**

**Results:**

|  |  |  |  |
|--|--|--|--|
| <b>Final burette reading(cm<sup>3</sup>)</b>   |  |  |  |
| <b>Initial burette reading(cm<sup>3</sup>)</b> |  |  |  |
| <b>Volume of FA3 used (cm<sup>3</sup>)</b>     |  |  |  |

(4  $\frac{1}{2}$  marks)

Record the titre values used to calculate average volume of **FA3** used.....cm<sup>3</sup>

( $\frac{1}{2}$  mark)

Average volume of **FA3**

used.....cm<sup>3</sup>  
(2  $\frac{1}{2}$  marks)

**Questions:**

(h) Calculate the number of moles;

(i) excess acid that reacted with hydroxide ions in **FA3**.

(1  $\frac{1}{2}$  marks)

.....

.....

.....

.....

.....

.....

(ii) excess acid in 250cm<sup>3</sup> of **FA5**.

(01 mark)

.....

.....

.....

(iii) acid that reacted with **Z**.

(02 marks)

.....

.....

.....

.....

.....

(c) Determine the reaction ratio between  $H_2X$  and  $Z$ .  
*(Relative formula mass of  $Z = 84$ )*

.....

.....

.....

.....

.....

.....

2. You are provided with substance **X** which contains **three** cations and **one** anion. You are required to carry out the tests in Table 3 below to identify the anion and cations in **X**. Identify any gas(es) evolved. Record your observations and deductions in the table. (30 marks)

**Table 3**

|     | TEST   | OBSERVATIONS | DEDUCTIONS |
|-----|--|--------------|------------|
| (a) | Heat a spatula endful of <b>X</b> in a dry test tube until there is no further change.   |              |            |
| (b) | Shake <b>two</b> spatula endfuls of <b>X</b> in a boiling tube with about $3\text{cm}^3$ of water. Add dilute sodium hydroxide solution to the mixture dropwise until in excess. |              |            |

|       |  |  |  |
|-------|--|--|--|
|       | Warm and filter. Keep both the filtrate and residue.   |  |  |
| (c)   | To the filtrate, add dilute nitric acid dropwise until the solution is just acidic. Divide acidic filtrate into <b>six</b> portions. |  |  |
| (i)   | To the <b>first</b> part of the filtrate, add sodium hydroxide solution dropwise until in excess.                                    |  |  |
| (ii)  | To the <b>second</b> part of the acidified filtrate, add ammonia solution dropwise until in excess.                                  |  |  |
| (iii) | To the <b>third</b> part of the filtrate add 2-3 drops of dilute sulphuric acid.   |  |  |
| (iv)  | To the <b>fourth</b> part of the filtrate, add 2-3 drops of litmus solution followed by ammonia solution dropwise until in excess.   |  |  |
| (v)   | To the <b>fifth</b> part, add lead(II) nitrate solution and heat.  |  |  |

|       |  |  |  |
|-------|--|--|--|
| (vi)  | <p>Use the <b>sixth</b> part of the acidified filtrate to carry out a test of your own to confirm the anion in <b>X</b>.</p> <p><b>Test:</b></p> <p>.....</p> <p>.....</p> <p>.....</p> <p>.....</p> |  |  |
| (d)   | <p>Wash the residue from part (b) with water and dissolve it dilute hydrochloric acid and divide the solution into <b>four</b> parts.</p>  |  |  |
| (i)   | <p>To the <b>first</b> part add sodium hydroxide solution dropwise until in excess.</p>  |  |  |
| (ii)  | <p>To the <b>second</b> part add dilute ammonia solution dropwise until in excess.</p>   |  |  |
| (iii) | <p>To the <b>third</b> part add, 4 drops of dilute hydrochloric acid, the 2 cm of magnesium ribbon and boil. Leave mixture to stand for 5 minutes. Then add potassium</p>                            |  |  |



|      |  |  |  |
|------|--|--|--|
|      | hexacyanoferrate(III) to the resultant solution.   |  |  |
| (iv) | To the fourth part of the acidified filtrate, add 3-4 drops of potassium thiocyanate solution. |  |  |

(f) Identify the anion and cations in **X**

Anion; .....

Cations; ....., ..... and.....

3. You are provided with organic substance **Y**. You are required to identify the nature of **Y**. carry out the tests below and record your observations and deductions. (20 marks)

|     | TESTS  | OBSERVATIONS | DEDUCTIONS |
|-----|--|--------------|------------|
| (a) | Burn a small amount of <b>Y</b> in a porcelain dish or on a spatula end  |              |            |
| (b) | Shake 0.5 cm <sup>3</sup> of <b>Y</b> with about 1 cm <sup>3</sup> of water and test the resultant solution with litmus paper. |              |            |
| (c) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add neutral iron(III) chloride solution.   |              |            |
| (d) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add a few drops of bromine water and shake.  |              |            |
| (e) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add potassium dichromate   |              |            |

|     |  |  |  |
|-----|--|--|--|
|     | solution and heat.   |  |  |
| (f) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add Lucas reagent.   |  |  |
| (g) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add 2-3 drops of Brady's reagent   |  |  |
| (h) | To about 0.5cm <sup>3</sup> of <b>Y</b> , add 5 drops of concentrated sulphuric acid followed by a concentrated solution of sodium hydrogen sulphite and shake strongly.   |  |  |
| (i) | Dissolve about 0.5cm <sup>3</sup> of <b>Y</b> in about 1 cm <sup>3</sup> of methanol. To the resultant solution, add iodine solution followed by dilute sodium hydroxide solution dropwise until the solution is pale yellow. Warm the mixture and allow to stand. |  |  |

Comment on the nature of **Y**.

.....

.....

.....

.....

END

## CONFIDENTIAL

Each student is to be provided with;

Volumetric flask (250ml)

Burette

Pipette (only 25.0cm<sup>3</sup>)

2 conical flasks

8 test tubes

3 beakers (250ml each)

100 or 50ml measuring cylinder

Distilled water

100 cm<sup>3</sup> of **FA1**

100 cm<sup>3</sup> of **FA2**

100 cm<sup>3</sup> of **FA3**

3.0g of solid **Z**

3.0g of solid **X**

Phenolphthalein indicator

A spatula

1 filter paper

Litmus papers

**Source of heat**

**Weighing scale**

8 cm<sup>3</sup> of **Y**

**FA1** is 1M sulphuric acid solution.

**FA2** is 0.25M sodium hydroxide solution

**FA3** is 0.041M sodium hydroxide solution.

**Z** is sodium hydrogen carbonate

**X** is a mixture of  $Al_2(SO_4)_3$ ,  $Fe_2(SO_4)_3$  (Ferric sulphate) and  $(NH_4)_2SO_4$  in a ratio 3:2:1

**Y** is phenylmethylketone or phenylethanone.